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INSTRUCTIONS FOR USE

EQUIPMENT:

Technical Name: X-Ray Viewer Holder Kit

Brand: GNATUS

Trade Name: X-Ray Viewer Holder Kit (200x300mm)

Manufacturer/ Distributor:

GNATUS - EQUIPAMENTOS MÉDICO-ODONTOLÓGICOS LTDA.
Rod. Abrão Assed , Km 53+450m - Cx. Postal 782 CEP 14097-500
Ribeirão Preto - S.P. - Brasil
Fone +55 (16) 2102-5000 - Fax +55 (16) 2102-5001
C.N.P.J. 48.015.119/0001-64 - Insc. Est. 582.329.957.115
www.gnatus.com.br - gnatus@gnatus.com.br

Technical Duties: Gilberto Henrique Canesin Nomelini
CREA-SP: 0600891412

ATTENTION

For greater safety:

Read and understand all the instructions contained in these Instructions for Use before installing or operating this Equipment.

Note: These Instructions for Use must be read by all the operators of this Equipment.

INDEX

PRESENTATION OF MANUAL 13

IDENTIFICATION OF EQUIPMENT 15

MODULES, ACCESSORIES, OPTIONS AND MATERIALS OF CONSUMPTION 17

TECHNICAL SPECIFICATIONS 18

- Installation..... 18

- Dimension..... 19

CORRECTIVE AND PREVENTIVE MAINTENANCE AND PRESERVATION..... 22

- Cleaning and disinfection..... 22

WARRANTY OF EQUIPMENT 23

FINAL CONSIDERATIONS..... 23

DESCRIPTION OF THE EQUIPMENT

Dear Customer

Congratulations. You have made a good choice when you decided to buy a GNATUS QUALITY product comparable to the best products available in the World. This manual is a general presentation of your product and it will give you important details to help you to solve possible problems. Please, read it and keep this with you.

Identification

Technical Name: X-Ray Viewer Holder Kit

Brand: GNATUS

Trade Name: X-Ray Viewer Holder Kit (200x300mm)

DESCRIPTION OF THE EQUIPMENT

Description and characteristics of the equipment

X-Ray viewer holder especially designed to be installed in the column of the dental lights with Ø 51mm or Ø 38mm.

Has both horizontal and vertical movement, providing better ergonomics.

Structure made of steel, with surface treatment via nanotechnology. Highly shiny smooth paint with an epoxy base, polymerized in an oven at 250°C, resistant to corrosion and cleaning materials.

Indication of the equipment


This equipment is exclusively for dental use, having to be employed and handled by a capacitated person (professional duly regulated, as per the local legislation of the country) observing the instructions contained in this manual.

The user is obliged to only use the equipment in perfect conditions and protect himself/herself, patients and third parties against possible hazards.

Purpose of the equipment

The X-Ray Viewer Holder Kit has the purpose of fastening the x-ray viewer screen on the dental light column, providing practicality in the viewing of x-rays.


MODULES, ACCESSORIES, OPTIONS AND MATERIALS OF CONSUMPTION

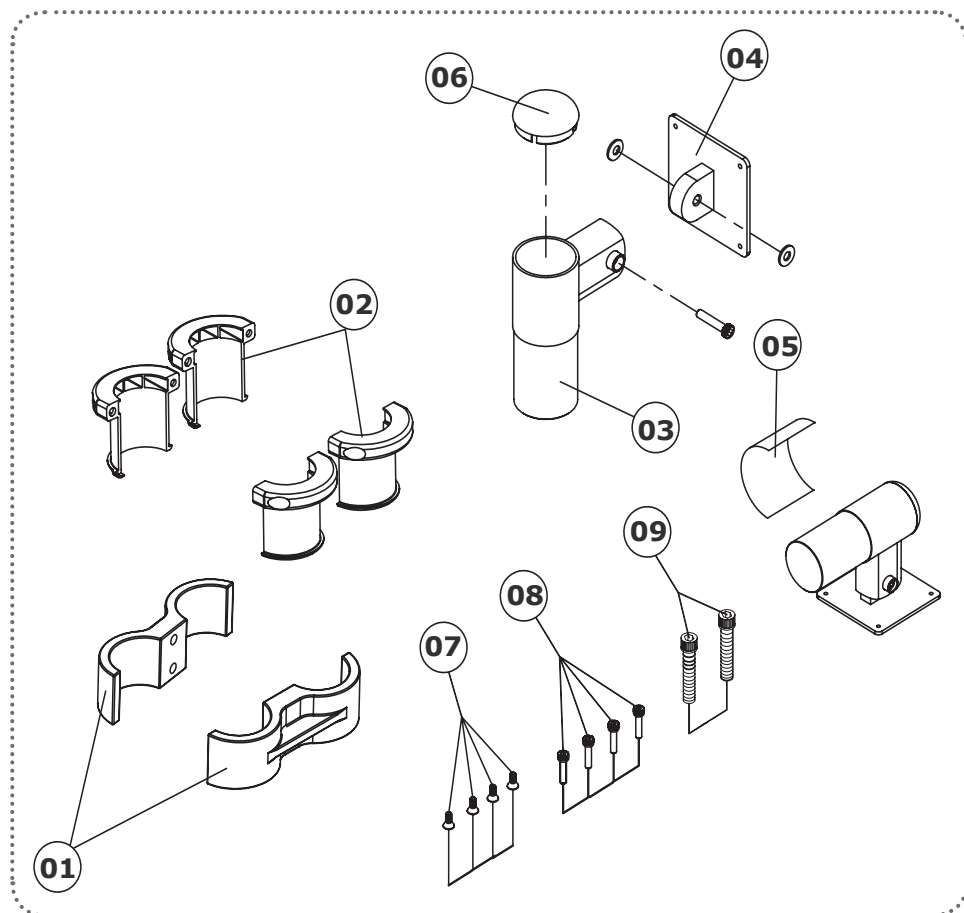- 01 - Unthreaded hanger adapter
- 02 - Internal brace
- 03 - S/C X-Ray Viewer Holder Articulation
- 04 - S/C X-Ray Viewer Holder Plate
- 05 - Courvin Vinyl Protector
- 06 - Finishing lid
- 07 - Steel Screw FL PH Mittopt 3x10mm head
- 08 - Allen Screw with M5X20mm head
- 09 - Allen Screw with M8x40MM head

TECHNICAL SPECIFICATIONS

Dimension (mm)

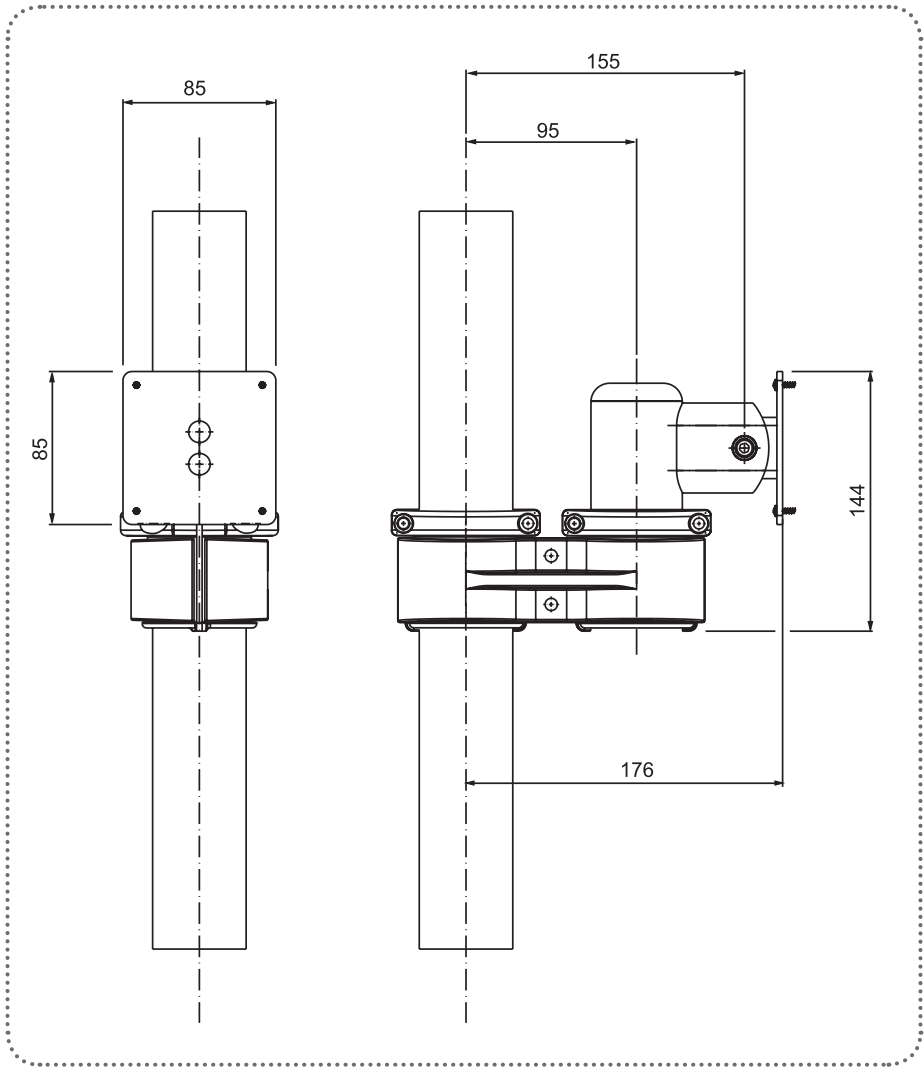
TECHNICAL SPECIFICATIONS

PLACE OF INSTALLATION OF THE X-RAY VIEWER HOLDER KIT


INSTALLATION

To fasten the x-ray viewer kit in the dental office, place the 4 braces (A) at the desired height around the dental light column. Then fit the 2 adaptors (B) over the braces and screw them on using the two Allen screws with M8 x 40mm head (C).

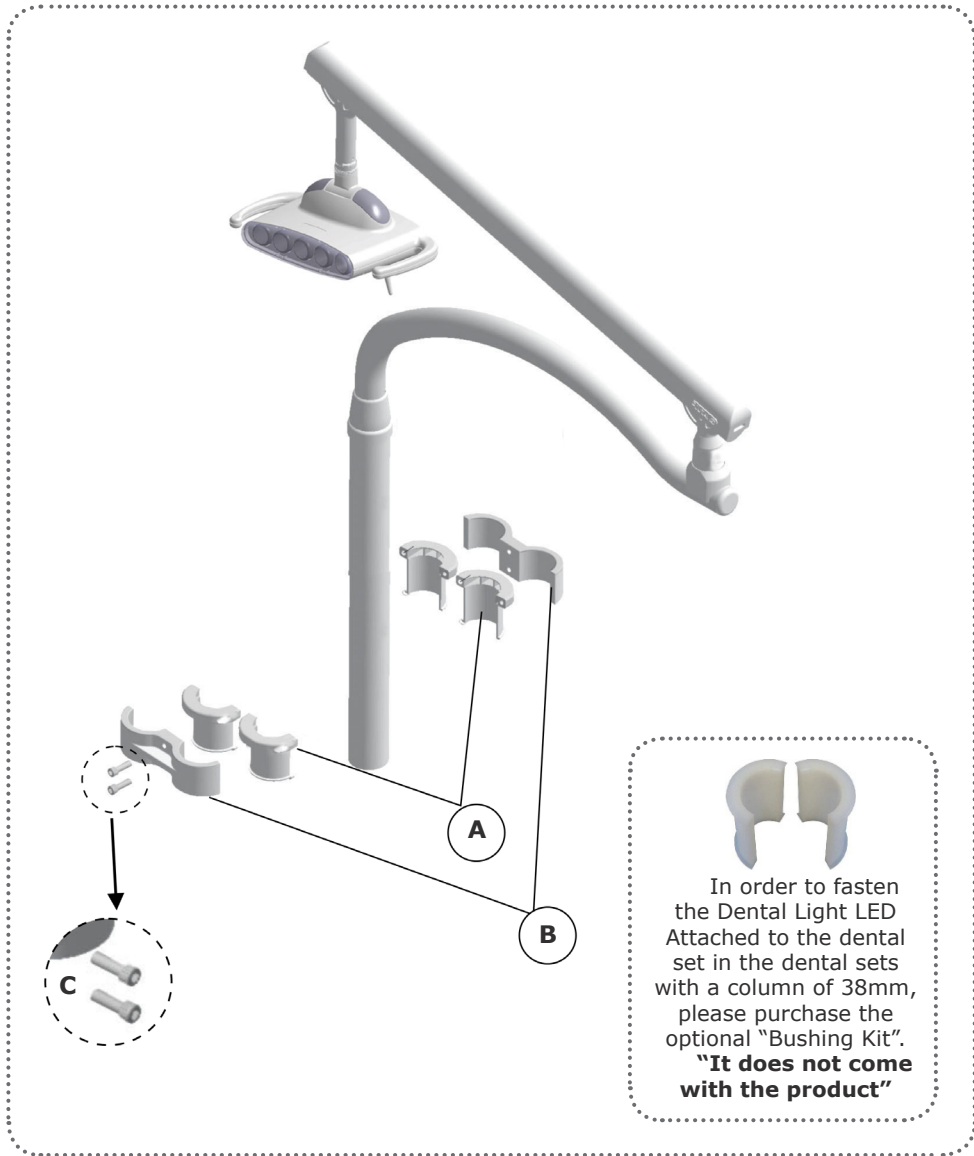INSTALLATION

Then fit the x-ray viewer holder (D) in the hangers (E).
Screw the x-ray viewer on the holder using the 4 steel screws FL PH 3x10mm head.

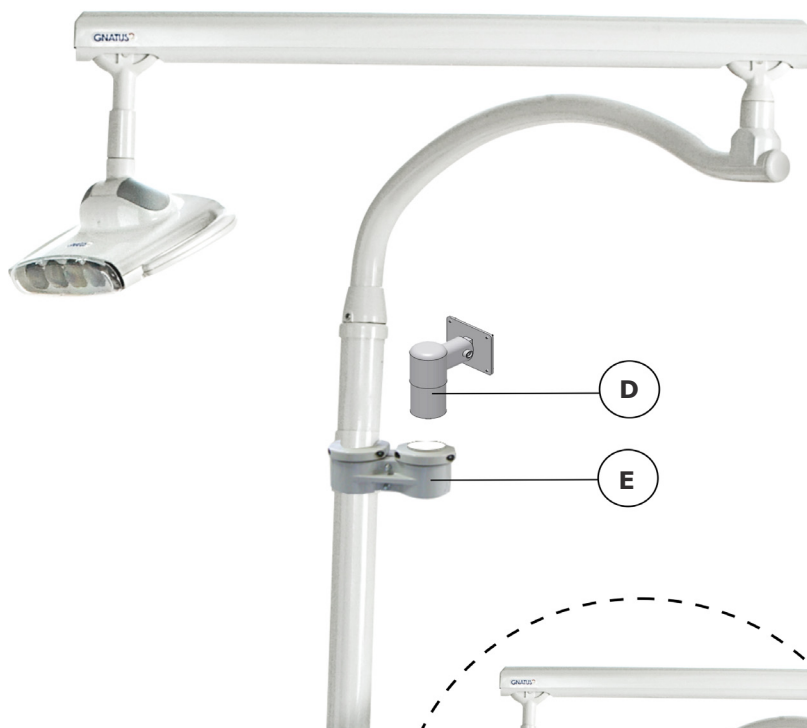CORRECTIVE AND PREVENTIVE MAINTENANCE AND PRESERVATION

Recommendations for the dental equipment maintenance.

Your Gnatus equipment has been designed and developed according to the standards of modern technology. Similarly to other kinds of equipment, it requires special care, which is many times neglected due to several reasons and circumstances.

Therefore, here are some important reminders for your daily routine. Try to follow these simple rules, which will save you a lot of time and will avoid unnecessary expenses once they start making part of your working procedure.

Additional procedures for reuse

The equipment can be reused in undetermined, i.e. unlimited, quantities, only needing to be cleaned and disinfected.

Cleaning

Chemical composition: Butyl Glycol, Decyl polyglucose, Sodium Benzoate, Sodium Nitrate, Essence, Deodorized Propane / Butane, demineralized Water.

For more information concerning cleaning procedures, see manufacturer's instructions.

WARNING:

- In order to prevent risks and damages to equipment, make sure that the liquid does not enter into the unit.
- The application of other solvent-based cleaning products or sodium hypochloride isn't recommended, because they may damage the equipment.

NOTE: The registration at the Ministry of Health of the "BactSpray" is executed separately from the product described in this manual, as the "BactSpray" is not manufactured by Gnatus.

Disinfection

Use clean and soft cloth dampened in alcohol 70% to disinfection of the equipment. Never use corrosive disinfectants or solvents.

Note: Use gloves and other systems of protection, during the disinfection.

Attention: Do not use any disinfectant spray, as the vapor may be inflammable, or it may cause injury.

WARRANTY OF EQUIPMENT

This equipment is covered by the warranty terms counting from the date of installation, as specified below; provided that the defect has occurred in normal conditions of use and that the equipment has not remained stored for more than 06 months counting from the issue date of the sales document until the date of the actual installation.

- WARRANTY TERMS: Verify the guarantee certificate;

- LOSS OF THE WARRANTY:

- A) Attempt to repair using an inadequate tool or by unauthorized technicians;

- B) Installation of the equipment by an unauthorized technician;

- C) Damage arising from inappropriate storage or signs of infringement;

- D) Incorrect use of the equipment;

- E) Use of a cleaning product not indicated by the factory;

- F) Falls or blows which the equipment may undergo or lack of observation of an compliance with the guidelines of the Owner's Manual, which was delivered with the present document, together with the equipment. Repair or replacement of parts during the warranty period shall not extend the validity term of their warranty.

- This warranty does not exempt the customer from paying the service charge for the visit and the travel expenses of the technician, except when the customer sends the equipment to execute the maintenance inside the establishment of the technical assistance.

- "Consumer Defense Code - art. 50, unique paragraph".

- The Warranty Certificate comes with the product and must be filled in upon the date of installation by the Gnatus Authorized Technician.

- Queries and information: GNATUS Help Desk (+55) 16 2102-5000.

- Check the warranty term attached to this manual.

ENGLISH

FINAL CONSIDERATIONS

The most important aspect related to equipment care is that concerning spare parts.

To guarantee the life span of your equipment, use only **original Gnatus spare parts**. They are sure to follow the technical specifications and standards required by Gnatus.

We must also point out to you our chain of authorized dealers. Only dealers that make part of this chain will be able to keep your equipment constantly new for they count on technical assistants who have been trained and on specific tools for the correct maintenance of your equipment.

Doubts and information: GNATUS Call center (55-16) 2102-5000.

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Obelis S.A, Boulevard Général Wahis 53, 1030 Brussels, Belgium,